

User Manual

1. Main page.....	2
2. Generic use for bacteria, yeast and primers.....	2
2.1. Entries.....	4
2.2. Filtering.....	5
2.2.1. Select User.....	5
2.2.2. Search for all fields.....	5
2.3. Ordering.....	6
2.4. Selecting different pages.....	7
2.5. Details of the data.....	8
2.6. Batch Create.....	9
2.7. Create/New.....	11
2.8. Edit.....	11
3. Data requirements.....	12
3.1. Bacteria.....	12
3.2. Primers.....	13
3.3. Yeast.....	13
4. Login and Sign out.....	14
4.1. Sign in.....	14
4.2. Sign out.....	14
5. Contact.....	14
6. Appendix.....	15
6.1. Appendix 1.....	15
6.2. Appendix 2.....	16
6.3. Appendix 3.....	16

1. Main page:

When you load in to the signed in main page, you will see a navigation column on the left-hand side of the page and pictures of the current lab in Manchester at the middle switching. (Figure 1)

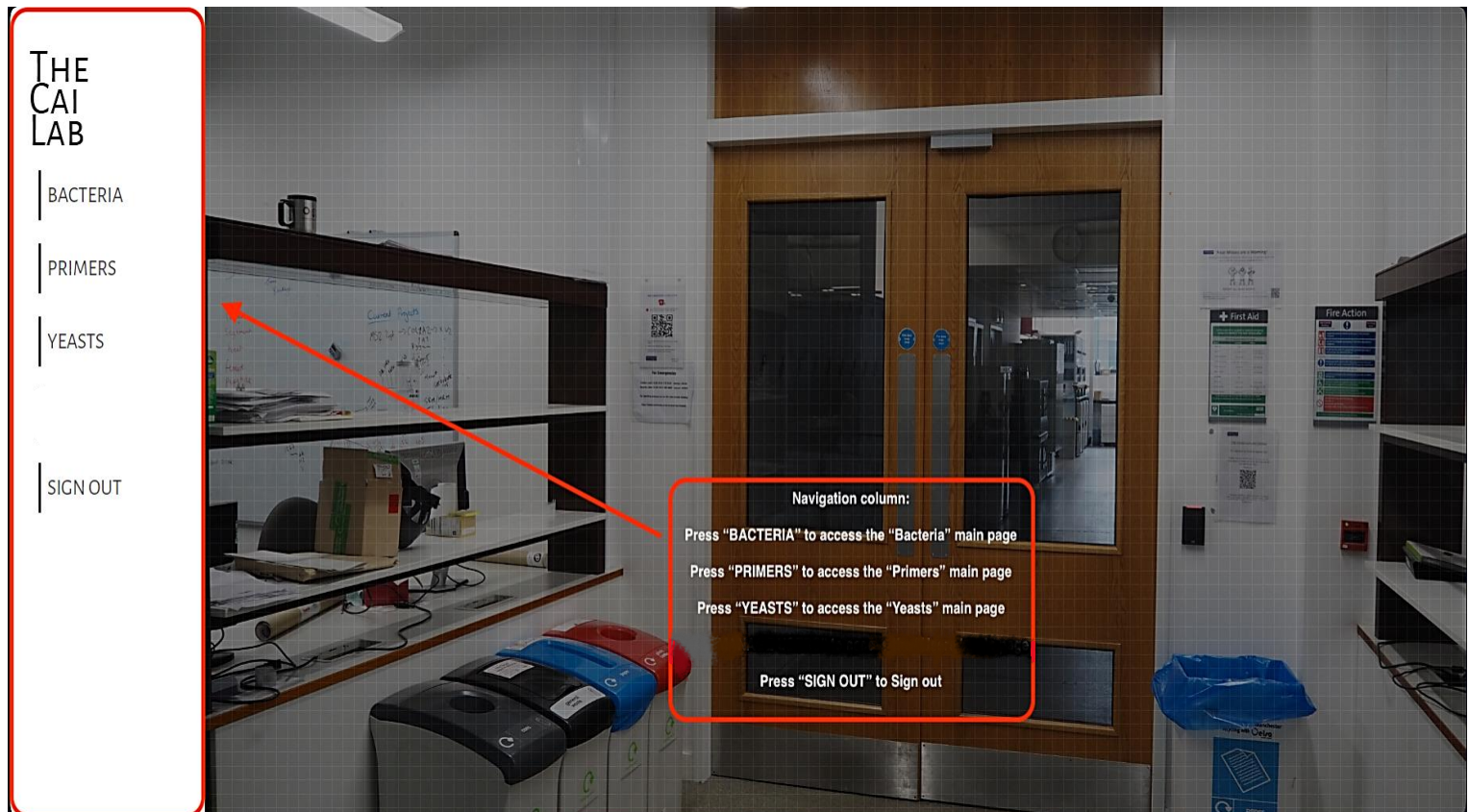


Figure 1: Signed in main page

2. Generic use for bacteria, yeast and primers:

The bacteria, primers, yeast pages follow the same structure. When you enter to either one of these pages, you will see a “select users” and “search all fields” filtering at the top left page under the title. You will see two buttons on the top right side, where it says, “Batch Create” and “Create New”. Under all of it, you will see a table where it lists out the data that exists in the database. (Figure 2, 3, 4)

Bacteria

Select Users:

Justin Leung (19)

Search for all fields in Bacteria:

Search bacteria

Batch Create

Create New Bacterium

Show 10 entries

Personal Id	Lab Id	Plasmid Name	Alternative Name	Host Strain	Comments	Location	Date
JLe019	YCe8010	pSCW11-Cre-EBD (URA3) test 17	pLM15, dae042	TOP10	FR18701117 from Jeff...	location B	2023-11-22 14:38:30
JLe018	YCe8009	pSCW11-Cre-EBD (URA3) test 16	pLM158, dae041	TOP10	FR18701117 from Jeff...	location B	2023-11-22 14:38:30
JLe017	YCe8008	pSCW11-Cre-EBD (URA3) test 17	pLM15, dae042	TOP10	FR18701117 from Jeff...	location B	2023-11-22 11:14:46
JLe016	YCe8007	pSCW11-Cre-EBD (URA3) test 16	pLM158, dae041	TOP10	FR18701117 from Jeff...	location B	2023-11-22 11:14:46
JLe015	YCe8006	pSCW11-Cre-EBD (URA3) test 17	pLM15, dae042	TOP10	FR18701117 from Jeff...	location B	2023-11-22 10:47:03
JLe014	YCe8005	pSCW11-Cre-EBD (URA3) test 16	pLM158, dae041	TOP10	FR18701117 from Jeff...	location B	2023-11-22

Figure 2: Bacteria Main page

Primers

Select Users:

Justin Leung (13)

Search for all fields in Primers:

Search primers

Batch Create

Create New Primer

Show 10 entries

Personal Id	Lab Id	Orientation	Melting Temperature	Concentration	Vendor	Location	Description
JLp013	YCp7188	Reverse	0				CDC27_40_target rv
JLp012	YCp7187	Forward	0				CDC27_40_target fw
JLp011	YCp7186	Reverse	12312312	12312312	12312312	12312312	12312312
JLp010	YCp7185	Forward	0	1	1	1	CDC27_40_target rv
JLp009	YCp7184	Forward	0	1	1	1	CDC27_40_target fw
JLp008	YCp7183	Forward	12312312	12312312	12312312	12312312	12312312
JLp007	YCp7182	Forward	0	1	1	1	CDC27_40_target rv
JLp006	YCp7181	Forward	0	1	1	1	CDC27_40_target fw
JLp005	YCp7180	Forward	0	1	1	1	CDC27_40_target rv

Figure 3: Primers Main page

Yeasts

Select Users:
Justin Leung (17) Search for all fields in Yeast:
Search yeast

Batch Create

Create New Yeast

Show 10 entries

Personal Id	Lab Id	Other Names	Mating Type	Plasmid	Genotype	Plasmid Type	Markers	Location	Comments	Parent Name	Other Parent Name	Date
JLy017	YCy5745	asdasd		asdasd	asdasd	asdasd	asdasd	asdasd	asdasd			2023-11-18 17:13:47.0
JLy016	YCy5744	BY4742	Alpha	pLeuOMeRS + pRS413-URA3_K209*	MAT Δ his3 Δ 1 leu2 Δ 0 lys2 Δ 0 ura3 Δ 0				Double Transformation...	JLy014	JLy014	2023-11-22 11:04:24.0
JLy015	YCy5743		Alpha	BY4742	MAT Δ his3 Δ 1 leu2 Δ 0 lys2 Δ 0 ura3 Δ 0		His3; Lys2; Ura3; Leu2		This is the yeast st...	JLy012	JLy014	2023-11-22 11:04:24.0
JLy014	YCy5742	BY4742		pLeuOMeRS + pRS413-URA3_K209*	MAT Δ his3 Δ 1 leu2 Δ 0 lys2 Δ 0 ura3 Δ 0				Double Transformation...			2023-11-22 10:57:20.0
JLy013	YCy5741		Alpha	BY4742	MAT Δ his3 Δ 1 leu2 Δ 0 lys2 Δ 0 ura3 Δ 0		His3; Lys2; Ura3; Leu2		This is the yeast st...	JLy012		2023-11-22 10:57:20.0
JLy012	YCy5740	BY4742		pLeuOMeRS	MAT Δ his3 Δ 1 leu2 Δ 0				Double			2023-11-

Figure 4: Yeast Main page

The following sections will give you detailed information regarding to the specific parts of the webpage. Note that most of the figures come from Bacteria pages.

2.1. Entries:

You can show how many entries to display with the select option located at the top left side above the table, you can choose 10, 25, 50 and 100 entries to display. The default choice is showing 10 entries per page. (Figure 5)

Bacteria

Select Users:
Justin Leung (19) Search for all fields in Bacteria:
Search bacteria

Batch Create

Create New Bacterium

Show 10 entries
10
25
50
100

Click on the select option between the text "Show" and "entries"

Choosing how many entries you want to display(10, 25, 50, 100)

Personal Id	Lab Id	Plasmid Name	Alternative Name	Host Strain	Comments	Location	Date
JLe019	YCe8010	pSCW11-Cre-EBD (URA3) test 17	pLM15, dae042	TOP10	FR18701117 from Jeff...	location B	2023-11-22 14:38:30
JLe018	YCe8009	pSCW11-Cre-EBD (URA3) test 16	pLM158, dae041	TOP10	FR18701117 from Jeff...	location B	2023-11-22 14:38:30
JLe017	YCe8008	pSCW11-Cre-EBD (URA3) test 17	pLM15, dae042	TOP10	FR18701117 from Jeff...	location B	2023-11-22 11:14:46
JLe016	YCe8007	pSCW11-Cre-EBD (URA3) test 16	pLM158, dae041	TOP10	FR18701117 from Jeff...	location B	2023-11-22 11:14:46

Figure 5: Entries function in Bacteria

2.2. Filtering:

There are two ways for filtering the table, you can select the user you want to see, and you can search all fields in the table. This is located between the show entries select option and title.

The following section will give you information how these two-filtering work.

2.2.1. Select user:

The options you can pick are all users and specific users, the number beside each specific users are the numbers of data that the specific user has created. The default selected user will be yourself, if you have created data within the specific table, if not the default selected user will be all users and it will not display your user within the selection box. (Figure 6)

THE CAI LAB

BACTERIA

PRIMERS

YEASTS

SIGN OUT

Bacteria

Select Users:

- Justin Leung (19)
- Justin Leung (19)
- All Users
- Aimee H (25)
- Alba Ballerini (296)
- Alba Ubide (156)
- Andrea Martella (555)
- Claire Chung (1)
- dabram a (44)
- Daniel Schindler (309)
- Elisa BM (4)
- Ellie Payne (21)
- Eloise O'Connor (129)
- Eva Garcia-Ruiz (982)
- Gabrielle David (5)
- Hugo Villanueva (17)
- James Sanders (110)
- Jamie Auxillos (1779)
- Joshua James (44)
- Junyang Ong (2)
- Katarina Jönsson (9)

Search for all fields in Bacteria:

Search bacteria

Batch Create

Create New Bacterium

Accession	Alternative Name	Host	Strain	Comments	Location	Date
CW11-11-22	pLM15, dae042	TOP10	FR18701117	location B	2023-11-22 14:38:30	
CW11-11-22	pLM158, dae041	TOP10	FR18701117	location B	2023-11-22 14:38:30	
CW11-11-22	pLM15, dae042	TOP10	FR18701117	location B	2023-11-22 11:14:46	

Click the select option below the "Select Users:" in the filtering section to display which user you want to display

Can pick yourself (If you have created data for this specific table), all users and other users

Figure 6: Select user in Bacteria

2.2.2. Search for all fields:

This field allows you to search anything from the table, it will search all data from the columns. You do not have to press the return/enter button to use this filter, the table will show words that you have typed in the search box automatically like a hot reload every time you have typed in a word. (Figure 7)

Bacteria

Select Users:

Justin Leung (19)

Search for all fields in Bacteria:

Search bacteria

Show 10 entries

Personal Id	Lab Id	Plasmid Name	Alternative Name	Host Strain	Comments	Location	Date
JLe019	YCe8010	pSCW11-Cre-EBD (URA3) test 17	pLM15, dae042	TOP10	FR18701117 from Jeff...	location B	2023-11-22 14:38:30
JLe018	YCe8009	pSCW11-Cre-EBD (URA3) test 16	pLM158, dae041	TOP10	FR18701117 from Jeff...	location B	2023-11-22 14:38:30
JLe017	YCe8008	pSCW11-Cre-EBD (URA3) test 17	pLM15, dae042	TOP10	FR18701117 from Jeff...	location B	2023-11-22 11:14:46

Just fill in the "Search for all fields in Bacteria" to search all columns in the table

No need to press return, it will auto-load every word you input

Batch Create

Create New Bacterium

Figure 7: Search all fields in Bacteria

2.3. Ordering:

The default order when you load the page are the dates, from the newest to the oldest. On each of the column title, you will see a grey diamond shaped button beside the title. (Figure 8) When you click on it the first time, it will order from the lowest to the highest or first alphabet to last alphabet from the specific column you have pressed from. (Figure 9) If you click it for the second time, the order of the specific column will be opposite. (Figure 10) However, you cannot do ordering with two columns at the same time.

Bacteria

Select Users:

Justin Leung (19)

Search for all fields in Bacteria:

Search bacteria

Show 10 entries

Personal Id	Lab Id	Plasmid Name	Alternative Name	Host Strain	Comments	Location	Date
JLe019	YCe8010	pSCW11-Cre-EBD (URA3) test 17	pLM15, dae042	TOP10	FR18701117 from Jeff...	location B	2023-11-22 14:38:30
JLe018	YCe8009	pSCW11-Cre-EBD (URA3) test 16	pLM158, dae041	TOP10	FR18701117 from Jeff...	location B	2023-11-22 14:38:30
JLe017	YCe8008	pSCW11-Cre-EBD (URA3) test 17	pLM15, dae042	TOP10	FR18701117 from Jeff...	location B	2023-11-22 11:14:46

Click the "diamond" shaped button to change the order of the table

Batch Create

Create New Bacterium

Figure 8: Ordering in Bacteria

Bacteria

Select Users: Search for all fields in Bacteria:

Show entries

Press the "diamond" button once to order the column from oldest to latest or first alphabet to last alphabet or first number to the last number

Personal Id	Lab Id	Plasmid Name	Alternative Name	Host Strain	Comments	Location	Date
ABe001	YCe6284	HcKan_P	HcKan_P	ccdB survival	HcKan_P vector used ...	Box1	2020-04-03 00:00:00
ABe002	YCe6285	HcKan_O	HcKan_O	ccdB survival	HcKan_O vector used ...	Box1	2020-04-03 00:00:00
ABe003	YCe6286	HcKan_T	HcKan_T	ccdB survival	HcKan_T vector used ...	Box1	2020-04-03 00:00:00
ABe004	YCe6287	POT1		ccdB survival	POT vector used in Y...	Box1	2020-04-03 00:00:00

Figure 9: Ordering in Bacteria when you click once

Bacteria

Select Users: Search for all fields in Bacteria:

Show entries

Press the "diamond" button twice to order the column from latest to oldest or last alphabet to first alphabet or last number to the first number

Personal Id	Lab Id	Plasmid Name	Alternative Name	Host Strain	Comments	Location	Date
YSe301	YCe5242	pRS415-YGR254W-YGR257C-YGR258C		DH5alpha	YGR254W, YGR257C and...		2018-03-19 00:00:00
YSe300	YCe4761	pCRCT-YGR139W		Beta-10	gRNA construct of YG...		2018-01-10 00:00:00
YSe299	YCe4760	pCRCT-YGR107W		Beta-10	gRNA construct of YG...		2018-01-10 00:00:00
YSe298	YCe4759	pCRCT-YGR076C		Beta-10	gRNA construct of YG...		2018-01-10 00:00:00

Figure 10: Ordering in Bacteria when you click twice

2.4. Selecting different pages:

At the bottom right of the table, you can access the next or previous page of the table if the table is not showing all the entries. It shows you the current page number and the next four-page number and the last page number. You can access the specific displayed page number by pressing on it. (Figure 11)

You can click the numbers to switch pages

Highlighted with blue box is your current page

You can click "Previous" or "Next" to also switch pages

		Cre-EBD (URA3) test 17			from Jeff...		11-22 10:47:03
JLe014	YCe8005	pSCW11-Cre-EBD (URA3) test 16	pLM158, dae041	TOP10	FR18701117 from Jeff...	location B	2023-11-22 10:47:03
JLe013	YCe8004	pSCW11-Cre-EBD (URA3) test 17	pLM15, dae042	TOP10	FR18701117 from Jeff...	location B	2023-11-22 09:39:45
JLe012	YCe8003	pSCW11-Cre-EBD (URA3) test 16	pLM158, dae041	TOP10	FR18701117 from Jeff...	location B	2023-11-22 09:39:45
JLe011	YCe8002	aaa	aaaaaa	aaa	aaaa	aaa	2023-11-10 14:40:49
JLe010	YCe8001	pSCW11-Cre-EBD (URA3) test 16	pLM158, dae041	TOP10	FR18701117 from Jeff...	location B	2023-11-17 14:09:30

Showing 1 to 10 of 7,823 entries

Previous 1 2 3 4 5 ... 783 Next

Figure 11: Selecting different pages

2.5. Details of the data:

To access the detailed information of a specific “Personal Id” or “Lab Id” on the table, you can simply click either the “Personal Id” or “Lab Id” number which is coloured blue with a line under it. (Figure 12)

Inside the page, you will see the “Lab Id” and “Personal Id” at the top of the page. Below them, you will see the detailed data from the specific “Lab Id” and “Personal Id”. In addition, there is also two buttons located at the top right, it is the “New” (“Create”) and “Edit” button. (Figure 13)

THE
CAI
LAB

BACTERIA

PRIMERS

YEASTS

SIGN OUT

Bacteria

Select Users:

Justin Leung (19)

Search for all fields in Bacteria:

Search bacteria

Show 10 entries

Personal Id	Lab Id	Plasmid Name	Alternative Name	Host Strain	Comments	Location	Date
JLe019	YCe8010	pSCW11-Cre-EBD (URA3) test 17	pLM15, dae042	TOP10	FR18701117 from Jeff...	location B	2023-11-22 14:38:30
JLe018	YCe8009	pSCW11-Cre-EBD (URA3) test 16	pLM158, dae041	TOP10	FR18701117 from Jeff...	location B	2023-11-22 14:38:30
JLe017	YCe8008	pSCW11-Cre-EBD (URA3) test 17	pLM15, dae042	TOP10	FR18701117 from Jeff...	location B	2023-11-22 11:14:46

Click the “Personal Id” or “Lab Id”(blue words with line under it) to access the detailed information for that particular “Personal Id” or “Lab Id”

Batch Create

Create New Bacterium

Figure 12: How to access the details of a bacteria via “Personal Id” and “Lab Id”

THE
CAI
LAB

BACTERIA
PRIMERS
YEASTS
SIGN OUT

JLe019, YCe8010

Lab Id	YCe8010
Personal Id	JLe019
Plasmid Name	pSCW11-Cre-EBD (URA3) test 17
Alternate Name	pLM15, dae042
Host Strain	TOP10
Comments	FR18701117 from Jeff lab Cre-EBD plasmid with URA3 marker
Location	location B
Marker	
Date	2023-11-22 14:38:30

Click “Edit” to edit this set of data from this page

Click “New” to create a new data for this specific table(In this picture, “Bacteria”)

New

Edit

Figure 13: Details of a bacteria page

2.6. Batch Create:

The “Batch Create” button is located at the top right of the webpage, it is beside the “Create New” button. (Figure 14) When you enter the “batch create” page, you will see a “choose file” button, a “Submit” button and “Download Batch Create Format” button. (Figure 15)

The “Download Batch Create Format” will give you a “csv” file with the structure to input the data for bacteria, primers or yeasts (Appendix 1, 2, 3). It is highly recommended for you to follow this structure and this file format for clean data in our database, or else the “batch create” will not accept your file and ask you to resubmit until you satisfy the requirements of the data requirements (Table 3.1, 3.2 , 3.3) and file format.

When you have finished inputting the data into the file, you may upload your file by clicking the “choose file” button, after that you click the “Submit” button to apply the batch create function. If your “batch create” file does not mee the data requirements (Table 3.1, 3.2, 3.3) or file requirements, it will pop up an error message stating which requirement you did not meet and ask you to resubmit it again until the requirements are met.

Bacteria

Select Users:

Justin Leung (19)

Search for all fields in Bacteria:

Search bacteria

Show 10 entries

Personal Id	Lab Id	Plasmid Name	Alternative Name	Host Strain	Comments	Location	Date
JLe019	YCe8010	pSCW11- Cre-EBD (URA3) test 17	pLM15, dae042	TOP10	FR18701117 from Jeff...	location B	2023- 11-22 14:38:30
JLe018	YCe8009	pSCW11- Cre-EBD (URA3) test 16	pLM158, dae041	TOP10	FR18701117 from Jeff...	location B	2023- 11-22 14:38:30
JLe017	YCe8008	pSCW11- Cre-EBD (URA3) test 17	pLM15, dae042	TOP10	FR18701117 from Jeff...	location B	2023- 11-22 11:14:46

Click "Create New Bacterium" to access the "Create" page for "Bacteria"

Click "Batch Create" to access the "Batch Create" page

Batch Create

Create New Bacterium

Figure 14: How to access Batch create and Create from main Bacteria page

Batch Create Bacterium

To upload Batch create file to do batch create

Choose file No file chosen

4. Click the "Submit" button to upload the file, if it does not meet the data requirements, it will bounce you back to this page with a pop up message specifying which data requirement(s) you did not meet. If successful, it will return you to the main "bacteria" page for the "bacteria" as an example

Submit +

3. Upload your file by clicking the "Choose file" button, remember the file format must be "CSV"

To download the format

Download Bacteria Batch Create Format

1. Click "Download the Bacteria Batch Create Format" button to download the "Batch Create" format for the the specific table, for the "bacteria" as an example

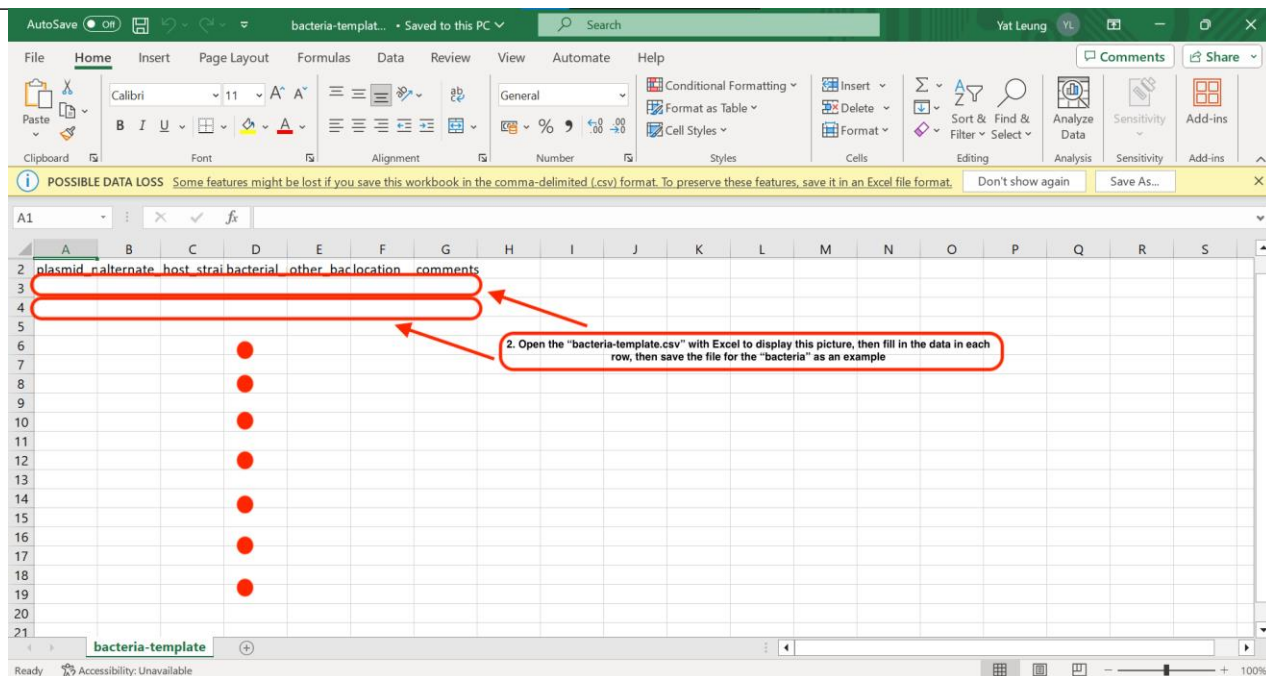


Figure 15: Batch Create Page for Bacteria

2.7. Create/New:

The “create” bacteria, primers and yeasts button is located at the top right side of the page beside the “Batch Create” button (Figure 14), or it is located top right in the details page of your specific “personal Id” and “lab Id” beside the “Edit” button, but is it labelled as “New”. (Figure 13)

Within the “create” page, you will see a form to ask you to fill in based on the specific data requirements (Table 3.1, 3.2, 3.3). After you have filled in the form, you can press the “Submit” button at the bottom right of the “create” page to create the data. However, if the data that you have input does not meet the requirements for the specific data requirement, it will pop up an error message stating which requirement(s) you did not meet on the “create” page, and ask you to resubmit it again until you meet the requirements. (Figure 16)

The screenshot shows the 'Create Bacterium' form on a web interface. On the left is a sidebar with 'THE CAI LAB' logo and navigation links for 'BACTERIA', 'PRIMERS', 'YEASTS', and 'SIGN OUT'. The main form is titled 'Create Bacterium' and contains several input fields: 'Plasmid name *', 'Alternate name', 'Comments' (a large text area), 'Date *', 'Host Strain', 'Location', 'Bacterial marker', and 'Other bacterial marker'. Below the form, there are two numbered instructions in red boxes: '1. Fill in the details for Bacteria based on the bacteria data requirements' and '2. Press “Submit” to create bacteria'. A red arrow points from instruction 1 to the form, and a green arrow points from instruction 2 to a 'Submit +' button.

Figure 16: Create Page for Bacteria

2.8. Edit:

The “edit” bacteria, primers and yeasts button is located top right in the details page of your specific “Personal Id” and “Lab Id” beside the “New” (“Create”) button. (Figure 13)

Within the “edit” page, you will see a form to ask you to edit based on the specific data requirements (Table 3.1, 3.2, 3.3). After you have edited the form, you can press the “Submit” button at the bottom right of the “edit” page to edit the data. However, if the data that you have input does not meet the requirements for the specific data requirement, it will pop up an error message stating which requirement(s) you did not meet on the “edit” page and ask you to resubmit it again until you meet the requirements. (Figure 17)

The screenshot shows the 'Edit Bacterium' form. A red box highlights the main form area, and a green box highlights the 'Submit +' button. A red arrow points from the text '1. Make the changes' to the form, and a green arrow points from the text '2. Click "Submit" to submit the changes' to the 'Submit +' button.

Edit Bacterium

Plasmid name *
pSCW11-Cre-EBD (URA3) test 17
Choose file No file chosen

Alternate name
pLM15, dae042

Comments
FR18701117 from Jeff lab Cre-EBD plasmid with URA3 marker

Date *
2023-11-22 14:38

Host Strain
TOP10

Location
location B

Bacterial marker

Other bacterial marker

1. Make the changes

2. Click "Submit" to submit the changes

Submit +

Figure 17: Edit page for Bacteria

3. Data requirements:

Table 3.1 Bacteria:

Name(UI Create and Edit Page)	Type	Null	Default	Comments
Plasmid Name	Text(255 characters)	No	None	
Alternate Name	Text(255 characters)	Yes	NULL	
Host Strain	Text(255 characters)	Yes	NULL	
Comments	Text	Yes	NULL	
Location	Text(255 characters)	Yes	NULL	
Date	Timestamp	No	NULL	
Bacterial marker	Int	Yes	NULL	Selective Option on Create and Edit page
Other bacterial marker	Text(255 characters)	Yes	NULL	
Choose file	Text(255 characters)	Yes	NULL	File type in .DNA

Table 3.2 Primers:

Name(UI Create and Edit Page)	Type	Null	Default	Comments
Description	Text	Yes	NULL	
Sequence	Text	Yes	NULL	
Melting temperature	Int	Yes	NULL	
Concentration	Text(255 characters)	Yes	NULL	
Vendor	Text(255 characters)	Yes	NULL	
Location	Text(255 characters)	Yes	NULL	
Date	Timestamp	No	CURRENT_TIMESTAMP	Default Generated, but must select a date on the edit and create page
Comments	Text	Yes	NULL	
Orientation	Int	Yes	NULL	Selective Option on Create and Edit page
Plate ID	Int	Yes	NULL	
Well ID	Text(255 characters)	Yes	NULL	

Table 3.3 Yeast:

Name(UI Create and Edit Page)	Type	Null	Default	Comments
Other names	Text(255 characters)	Yes	NULL	
Plasmid	Text(255 characters)	Yes	NULL	
Genotype	Text(255 characters)	Yes	NULL	
Plasmid type	Text(255 characters)	Yes	NULL	
Location	Text(255 characters)	Yes	NULL	
Comments	Text	Yes	NULL	
Date	Timestamp	No	CURRENT_TIMESTAMP	Default generated
Other mating type	Text(255 characters)	Yes	NULL	
Mating type	Int	Yes	NULL	Selective Option on Create and Edit page
Parent Yeast Id	Text(255 characters)	Yes	NULL	Parent name must exist within table, or else it will reject request for creating or editing
Other Parent Yeast Id	Text(255 characters)	Yes	NULL	Parent name must exist within table, or else it will reject request for creating or editing

Yeast Markers	Text(255 characters)	Yes	NULL	
---------------	----------------------	-----	------	--

4. Login and Sign out:

4.1. Sign in:

The “SIGN-IN” button is located at the left navigation column below the “WEBSITE” button on the before signed in main page. It will ask you to input your “Username” with your registered email and your “Password” that you have registered with this webpage. To login, you must press the “Log in” button, if one of your fields are incorrect, it will ask you to refill the “Username” and “Password” to log in again.(Figure 26)

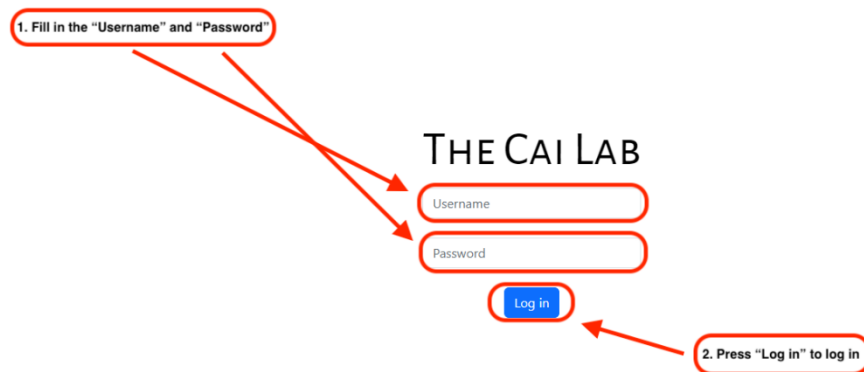


Figure 26: Sign in page

4.2. Sign out:

The “SIGN OUT” button is located at the left navigation column below “YEAST” button on the signed in pages. You simply need to click the “SIGN OUT” button to sign out.(Figure 1)


5. Contact:

If you have any queries or questions or run into any problems, please contact system administrator:

Raymond Wan: Raymond.wan@manchester.ac.uk


6.Appendix:

Appendix 1:

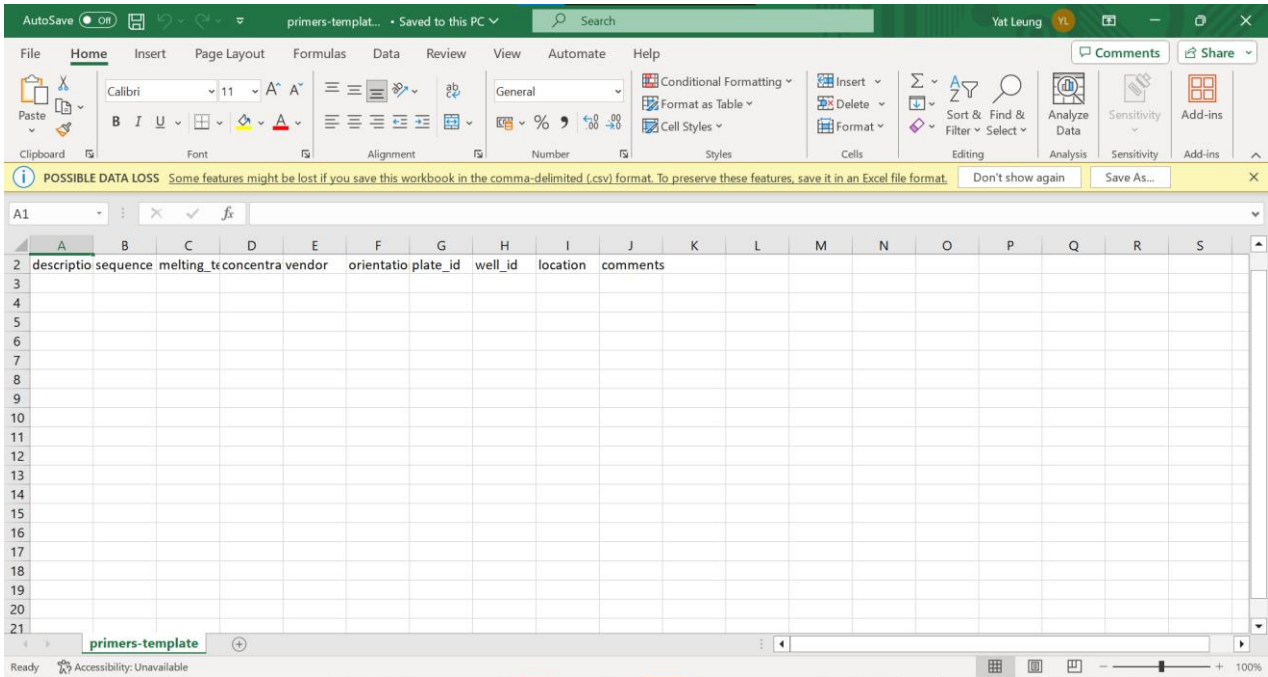


Appendix 1: Batch create structure for Bacteria

Appendix 2:

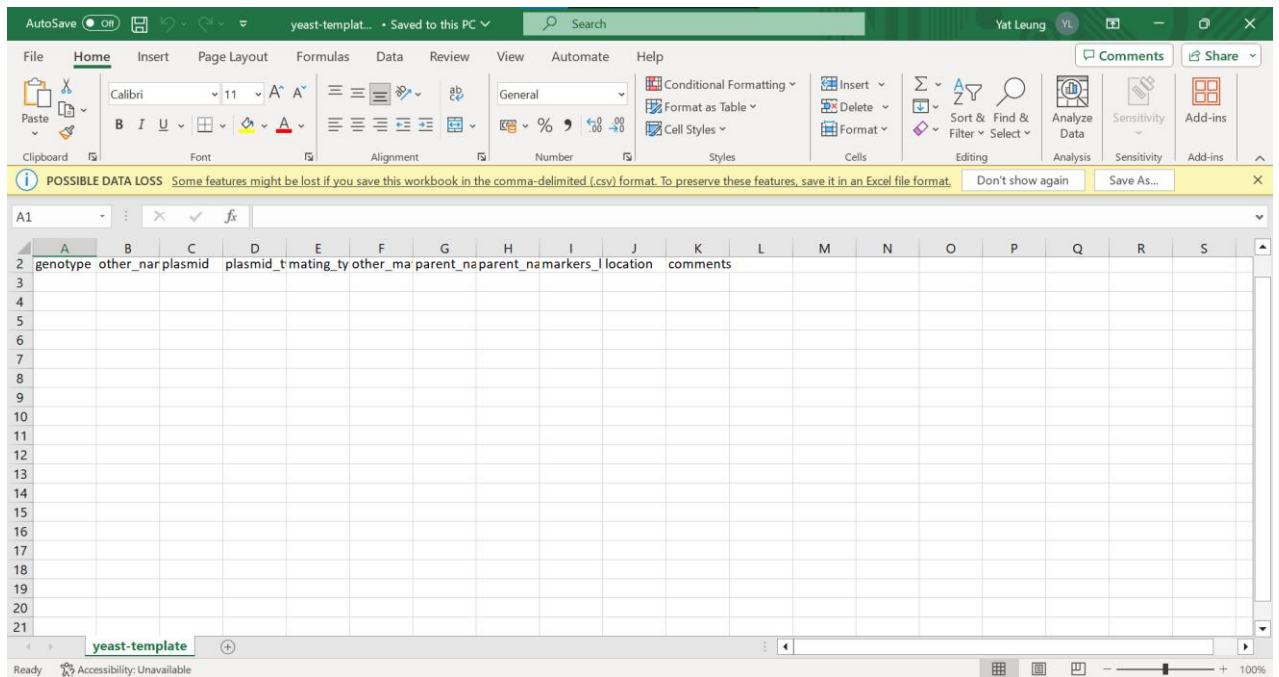

Appendix 2: Batch create structure for Primers

Appendix 3:

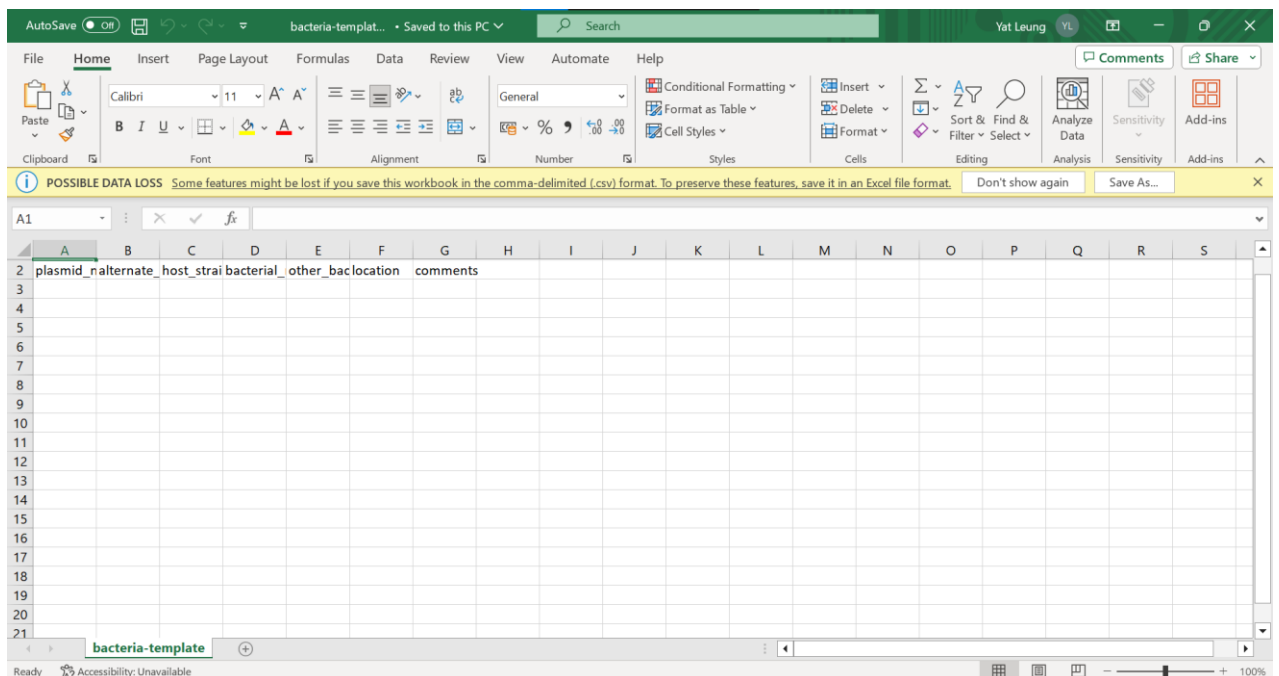

Appendix 3: Batch create structure for Yeast